





MOBILO FS

Functional Emergency Stretcher

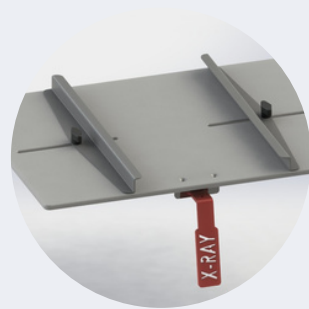

X-Ray Casette Holder

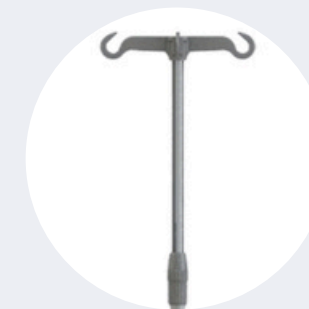

Height adjustable IV Pole Stand



Double Shock Absorber on Back




Emergency Shock Pedal



Oxygen Cylinder Holder



5th Castor

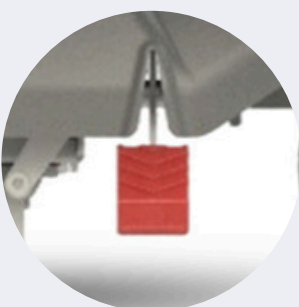

Central Locking Pedal



HPL Lying Surface



Urine Bag Hanger Hooks



Foldable Push Bar



Dış Ölçüler
External Dimensions
81 x 210 cm



Yatak Ölçüsü
Bed Dimensions
60,5 x 190,8 x 8 cm



Yükseklik (Max)
Height (Max)
90 cm



Yükseklik (Min)
Height (Min)
60 cm



Ters Trendelenburg
Reverse Trendelenburg
-15°



Trendelenburg
Trendelenburg
15°



Sırt Hareketi
Back Movement
90°



Teker Çapı
Castor Diameter
200 mm



Güvenli Çalışma Yüğü
Safe Working Load
320kg



MOBILO FS

Functional Emergency Stretcher

Genel Özellikler

- Nano seramik kaplı elektrostatik toz boyalı metal aksam
- Kompak malzemeden imal 4 parçalı yatma yüzeyi
- Kolay açılıp kapanabilen yan korkuluklar
- Dört köşesinde koruyucu tamponlar
- Ayarlanabilen katlanır serum askısı
- Plastik kaplı katlanabilen itme barları
- 360° dönebilen merkezi fren sistemli tekerler (200 mm)
- 5.Teker (Sedyede istenilen yönde hareket kolaylığı sağlar)
- Ayak pedalları ile hidrolik kontrollü yükseklik, trendelenburg ve ters trendelenburg ayarı
- Gazlı amartisör kontrollü sırt ayarı
- Diz hareketlerini kademeli mekanizma ile ayarlayabilme (Opsiyonel)
- ABS plastik kaplı alt şase
- Oksijen tüpü taşıma yeri
- 28 dansite şilte
- X-ray kaset taşıyıcısı
- Acil şok pozisyonu

General Features

- Nano ceramic coated electrostatic powder painted metal frame
- 4-piece lying surface made of compact material
- Side rails that can be easily opened and closed
- Protective bumpers on all four corners
- Adjustable folding IV pole
- Plastic-coated foldable push bars
- 360° swivel wheels with central brake system (200 mm)
- 5th wheel (allows easy movement in the desired direction on the stretcher)
- Hydraulic height, Trendelenburg, and reverse Trendelenburg adjustment with foot pedals
- Gas-filled shock absorber-controlled backrest adjustment
- Knee movement adjustment with a step mechanism (Optional)
- ABS plastic-coated lower frame
- Oxygen cylinder storage compartment
- 28 density mattress
- X-ray cassette holder
- Emergency shock position



90° Sırt Hareketi

90° Backrest Movement



60-90 cm Yükseklik Hareketi

60-90 cm Height Movement



15° Trendelenburg

15° Trendelenburg



-15° Ters Trendelenburg

-15° Reverse Trendelenburg